



Custom Solutions




FENESTRATED TREO[®]

Advancing Fenestrated Horizons


Trimodular System

The Fenestrated TREO platform, based on the TREO Abdominal Stent-Graft System, is a simple trimodular system.

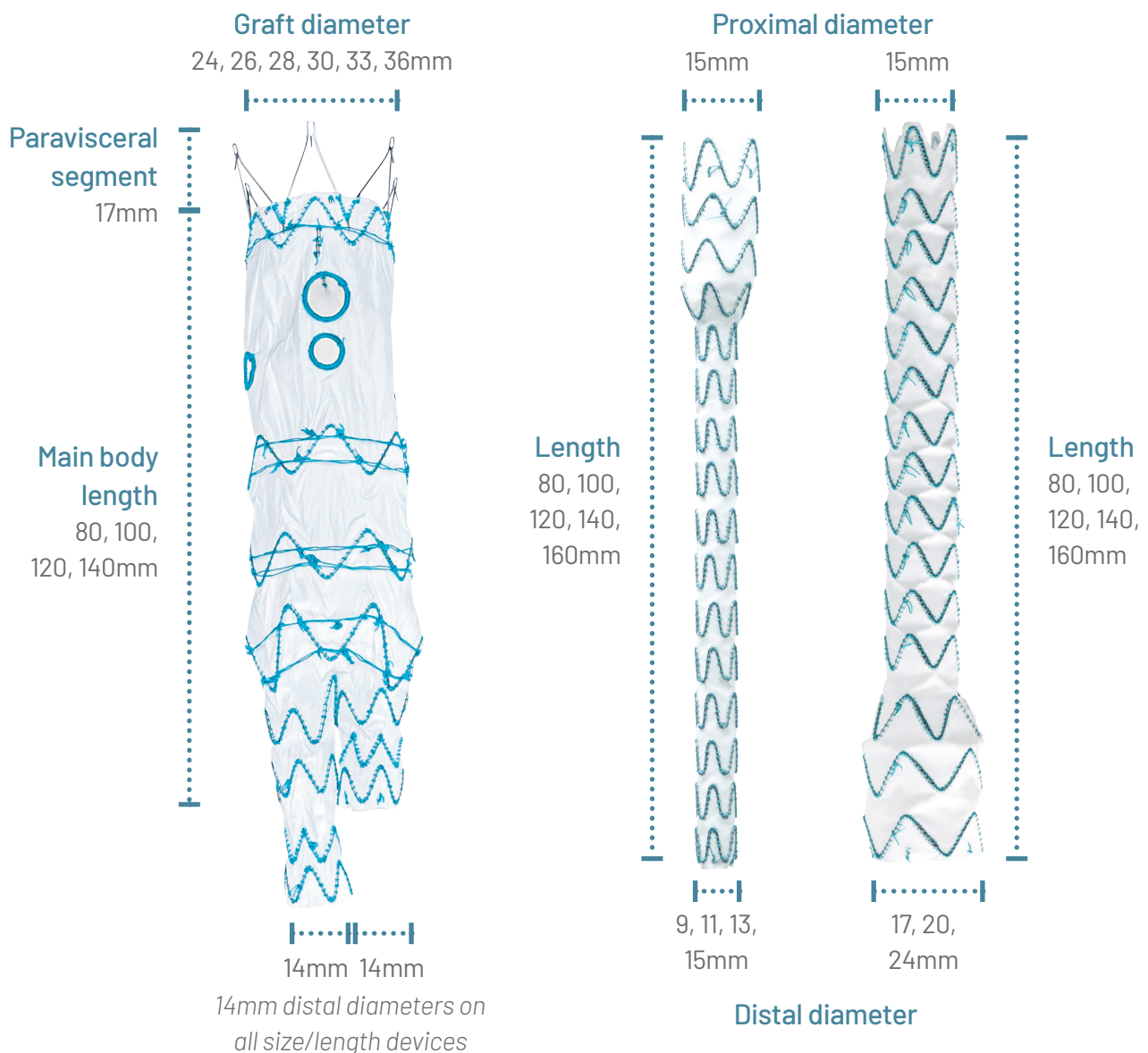

Reliable staged expansion for a controlled deployment*

Vessel diameter in the visceral area	Graft diameter	Constrained graft diameter
20mm	24mm	16mm
20mm	26mm	16mm
22mm	28mm	18mm
23mm	30mm	19mm
26mm	33mm	22mm
28mm	36mm	24mm

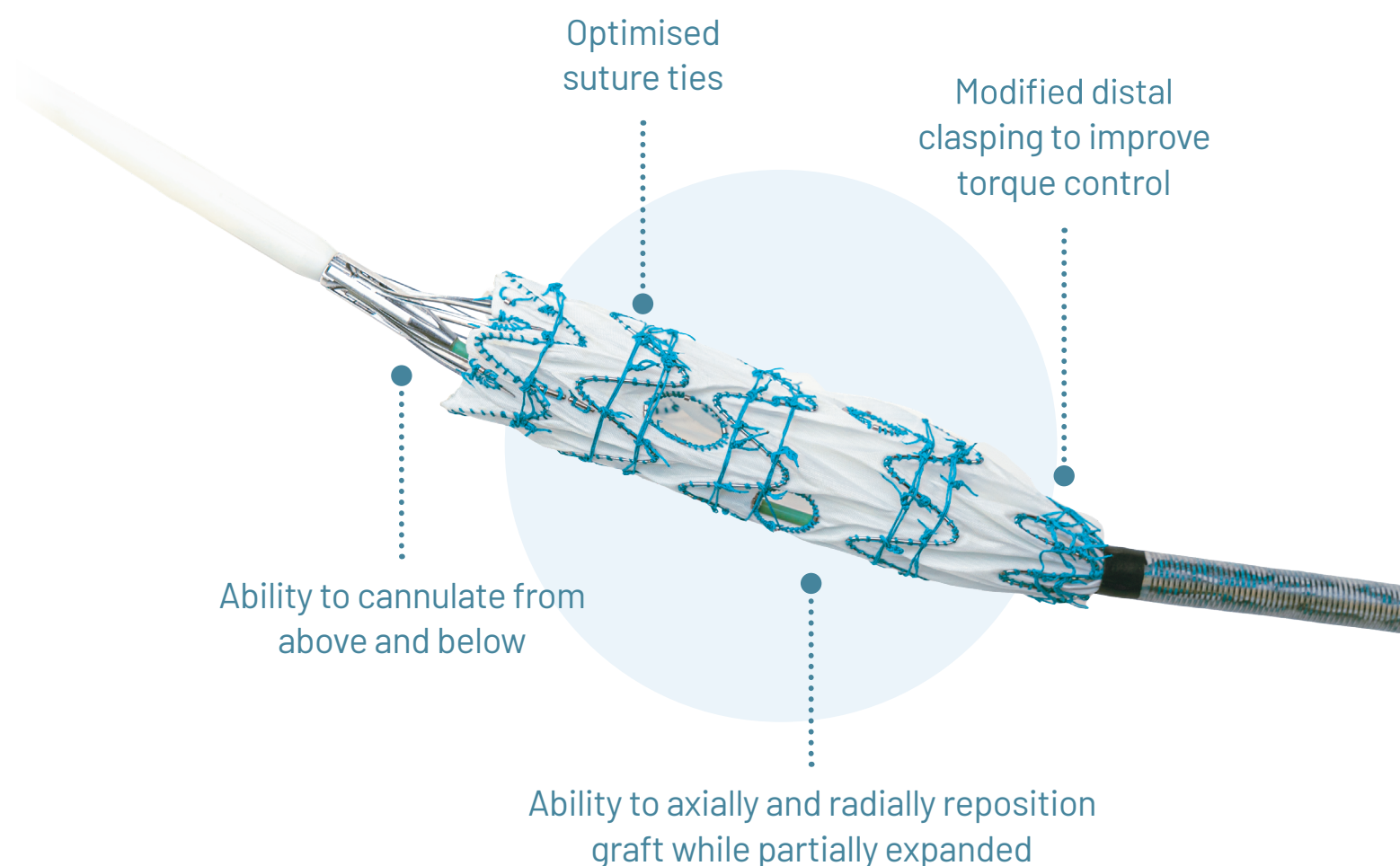

*Based on TREO Abdominal Stent-Graft System



Freedom to Position Fenestrations

- ▶ **Up to 5 fenestrations**
- ▶ Freedom to position stents and fenestrations as required
- ▶ **Reinforced fenestrations** with **Nitinol rings**
- ▶ **Circumferential radiopaque marker** for improved visualisation and alignment accuracy
- ▶ Main body can allow up to 40mm of free space between stents for fenestration placement
- ▶ **Fenestration sizes up to 12mm**



Optimised Stent Configuration

- ▶ Fully configurable support positioning for optimal fenestration placement
- ▶ Additional sealing stents can be added to provide a longer proximal sealing zone





Fully supported main body





Partially supported main body

Ability to Cannulate from Above

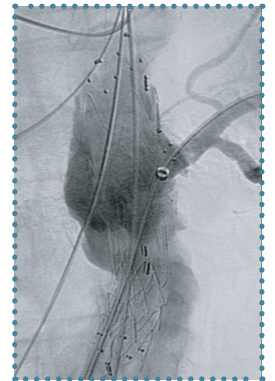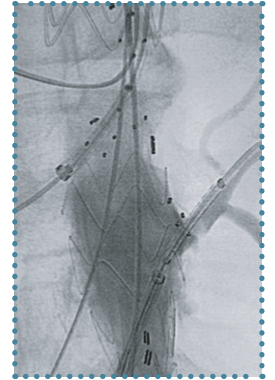

Circumferential radiopaque marker for
improved visualisation and alignment accuracy

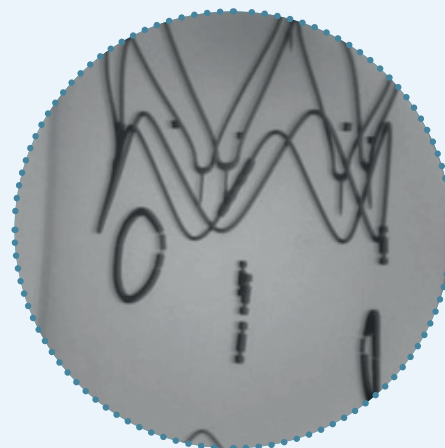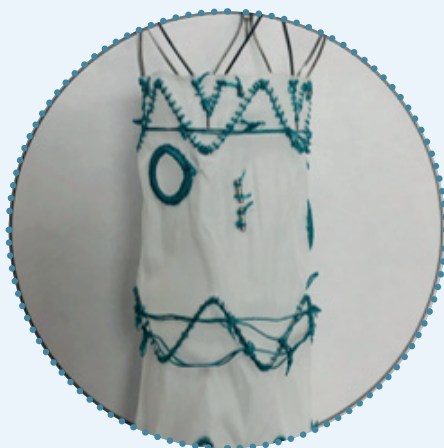

Trusted Low Profile Device*

19Fr

Main body

Leg extension

24 to 36mm

9 to 15mm

17 to 24mm

19 Fr(OD)

13 Fr(OD)

14 Fr(OD)

Reducing ties

Reducing ties mechanism allows graft repositioning until completely opened

Introducer sheath

60cm low profile sheath with hydrophilic coating and flexible tip for easier navigation


Precise delivery system

The mechanical deployment provides controlled and stable stent-graft deployment

Proximal clasp

The clasp mechanism keeps control of deployment and allows precise placement with cranial and caudal adjustment before the bare stent is released

Excellent torque control: distal leg clasp

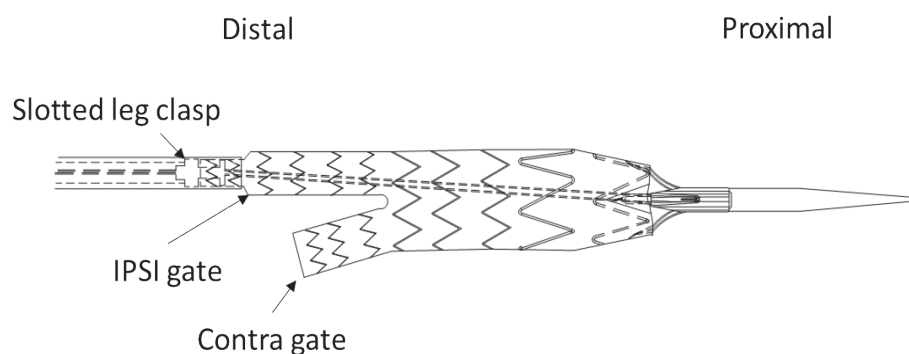

5 weeks
delivery
time^Δ

* Based on TREO Abdominal Stent-Graft System

^Δ Delivery time given as an estimate only. Please check current delivery times with your representative.





At Terumo Aortic, we are 100% focused on
addressing every segment of the aorta,
from the aortic root to the iliacs.

terumoaortic.com

Discover solutions for every segment of the aorta



Visit our website for more information on use, indications, contraindications, warnings/precautions and availability within your market.

Custom made devices are specifically made in accordance with a written prescription of any person authorised by national law by virtue of that person's professional qualifications; which gives (1) specific design characteristics provided under that person's responsibility and (2) is intended for the sole use of a particular patient exclusively to meet their individual conditions and needs. Custom made devices are not available in the US and availability is subject to local regulatory approval.

 Manufactured by: Bolton Medical Inc, 799 International Parkway, Sunrise, Florida 33325, United States

OUS_0249-A Mar23 • PM-06747

Distribuido por:

cardiva[®]

Fabricado por:

**TERUMO[®]**
Aortic